

Effective Date: Monday, May 02, 2016

Test Updates


In our continuing effort to provide you with the highest quality toxicology laboratory services available, we have compiled important changes regarding a number of tests we perform. Listed below are the types of changes that may be included in this notification, effective Monday, May 02, 2016

Test Changes -	Tests that have had changes to the method/ CPT code, units of measurement, scope of analysis,								
reference comments, or specimen requirements.									
Discontinued Te	sts - Tests being discontinued with alternate testing suggestions.								

Please use this information to update your computer systems/records. These changes are important to ensure standardization of our mutual laboratory databases.

If you have any questions about the information contained in this notification, please call our Client Support Department at (866) 522-2206. Thank you for your continued support of NMS Labs and your assistance in implementing these changes.

The CPT Codes provided in this document are based on AMA guidelines and are for informational purposes only. NMS Labs does not assume responsibility for billing errors due to reliance on the CPT Codes listed in this document.

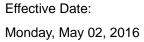



Test Updates

Test	Test Name	Test Name	TAT	Method / CPT Code		Stability	Analytes	Reference Comments	Discontinue
0556U	Benzidine, Urine				•	•			
1549U	Dichlorobenzidine, Urine				•	•			
3081U	MOCA and MDA Exposure Profile, Urine				•	•			
3080U	MOCA, Urine				•	•			
2986U	Methylenedianiline, Urine				•	•			



Effective Date:

Monday, May 02, 2016

Test Updates

Test Changes

0556U Benzidine, Urine

Summary of Changes: Specimen Requirements (Rejection Criteria) were changed.

Stability was changed.

Requested Volume: 5 mL Urine Transport Temperature: Frozen

Sample Container: Plastic container (preservative-free)

Light Protection: Not Required

Special Handling: None

Rejection Criteria: Received Room Temperature. Received Refrigerated.

Stability: Room Temperature: Not Stable

Refrigerated: 2 day(s) Frozen (-20 °C): 30 day(s)

1549U Dichlorobenzidine, Urine

Summary of Changes: Specimen Requirements (Transport Temperature) were changed.

Specimen Requirements (Sample Requirements) were changed. Specimen Requirements (Special Handling) were changed. Specimen Requirements (Rejection Criteria) were changed.

Stability was changed.

Requested Volume: 5 mL Urine
Transport Temperature: Refrigerated

Sample Container: Plastic container (preservative-free)

Light Protection: Not Required

Special Handling: None

Rejection Criteria: Received Room Temperature.

Stability: Room Temperature: 2 day(s)

Refrigerated: 30 day(s) Frozen (-20 °C): 30 day(s)

3081U MOCA and MDA Exposure Profile, Urine

Summary of Changes: Specimen Requirements (Requested Volume) were changed.

Specimen Requirements (Transport Temperature) were changed. Specimen Requirements (Sample Requirements) were changed. Specimen Requirements (Special Handling) were changed. Specimen Requirements (Rejection Criteria) were changed.

Stability was changed.



Effective Date:

Monday, May 02, 2016

Test Updates

Test Changes

Requested Volume: 5 mL Urine
Transport Temperature: Refrigerated

Sample Container: Plastic container (preservative-free)

Light Protection: Not Required

Special Handling: None

Rejection Criteria: Received Room Temperature.

Stability: Room Temperature: Not Stable

Refrigerated: 14 day(s) Frozen (-20 °C): 14 day(s)

3080U MOCA, Urine

Summary of Changes: Specimen Requirements (Special Handling) were changed.

Stability was changed.

Requested Volume: 5 mL Urine
Transport Temperature: Refrigerated

Sample Container: Plastic container (preservative-free)

Light Protection: Not Required

Special Handling: None

Rejection Criteria: Received Room Temperature.

Stability: Room Temperature: 2 day(s)

Refrigerated: 14 day(s) Frozen (-20 °C): 14 day(s)

2986U Methylenedianiline, Urine

Summary of Changes: Specimen Requirements (Transport Temperature) were changed.

Specimen Requirements (Sample Requirements) were changed. Specimen Requirements (Special Handling) were changed. Specimen Requirements (Rejection Criteria) were changed.

Stability was changed.

Requested Volume: 5 mL Urine
Transport Temperature: Refrigerated

Sample Container: Plastic container (preservative-free)

Light Protection: Not Required

Special Handling: None

Rejection Criteria: Received Room Temperature.

Stability: Room Temperature: Not Stable

Refrigerated: 14 day(s) Frozen (-20 °C): 14 day(s)